



ステロイド  
Steroids

X150801J

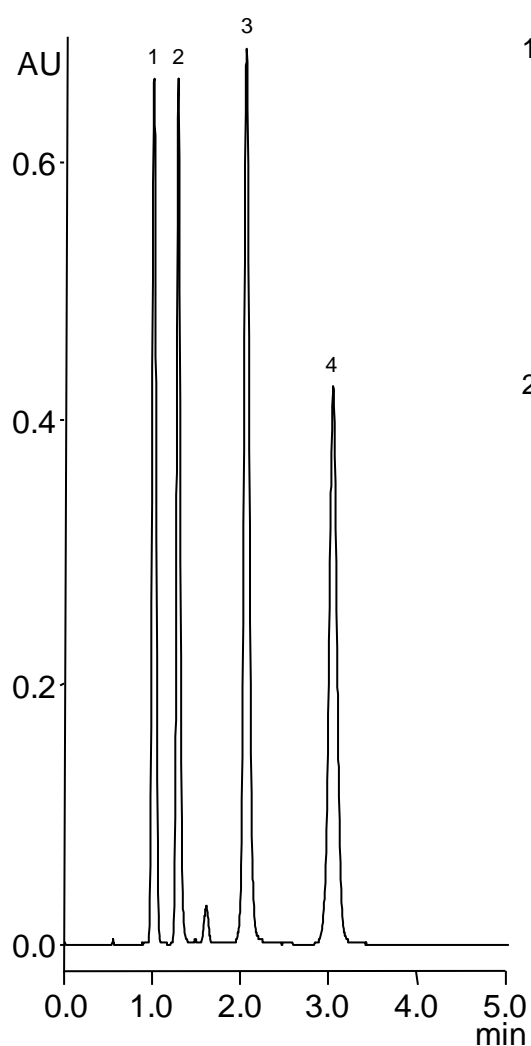


1.

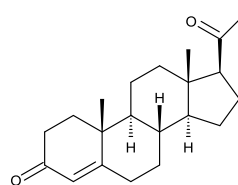


Progesterone

2.

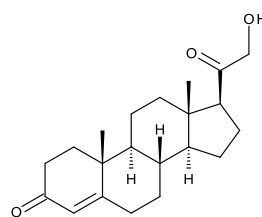


Deoxycorticosterone

3.

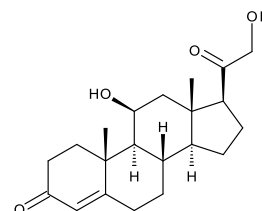

Corticosterone

4.

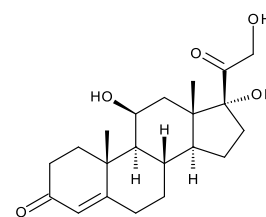

Hydrocortisone

Column : Alcyon SFC Triart Diol (5  $\mu$ m, 12 nm)  
150 X 4.6 mm I.D.  
Eluent : CO<sub>2</sub>/ethanol (80/20)  
Flow rate : 3.0 mL/min  
Temperature : 40°C  
Detection : UV at 254 nm  
Back Pressure : 13.8 MPa (2000 psi)  
Injection : 5  $\mu$ L (0.8 mg/mL)