



ステロイド
Steroids

X150803P

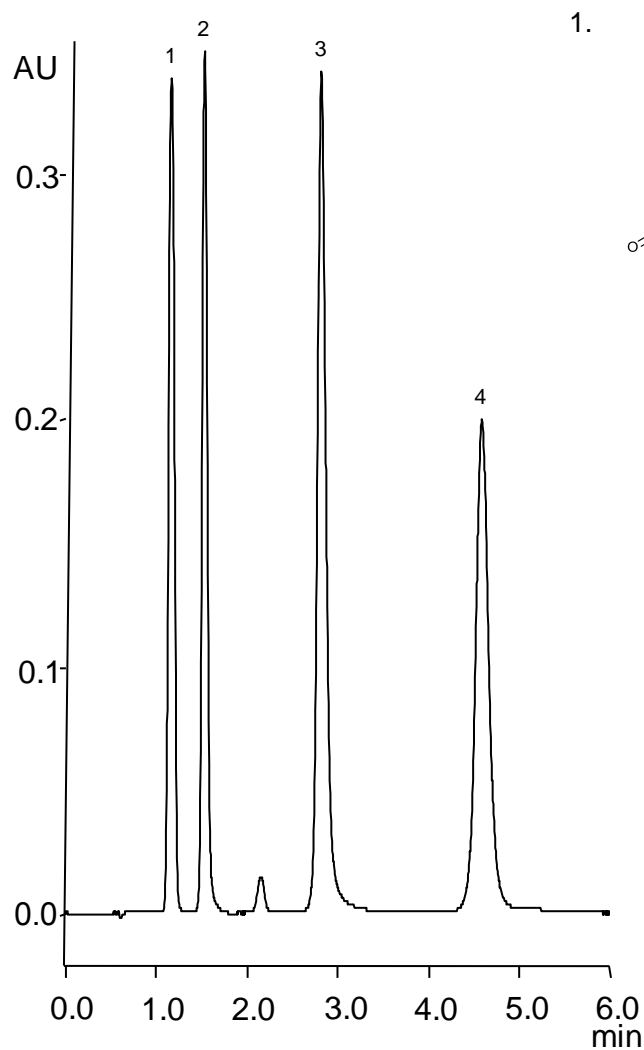


1.

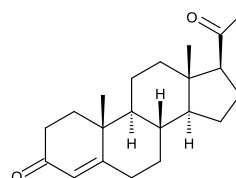


Progesterone

2.

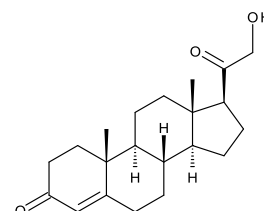


Deoxycorticosterone

3.

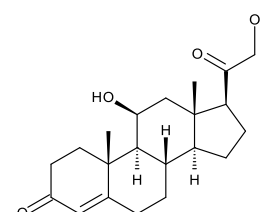

Corticosterone

4.

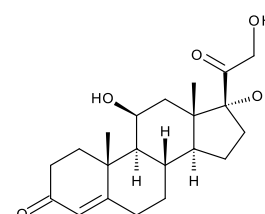

Hydrocortisone

Column	: Alcyon SFC Triart Diol (5 μ m, 12 nm) 150 X 4.6 mm I.D.
Eluent	: CO ₂ /2-propanol (80/20)
Flow rate	: 3.0 mL/min
Temperature	: 40°C
Detection	: UV at 254 nm
Back Pressure	: 13.8 MPa (2000 psi)
Injection	: 5 μ L (0.8 mg/mL)