



感冒藥

Ingredients in a cough/cold medication

X15117D

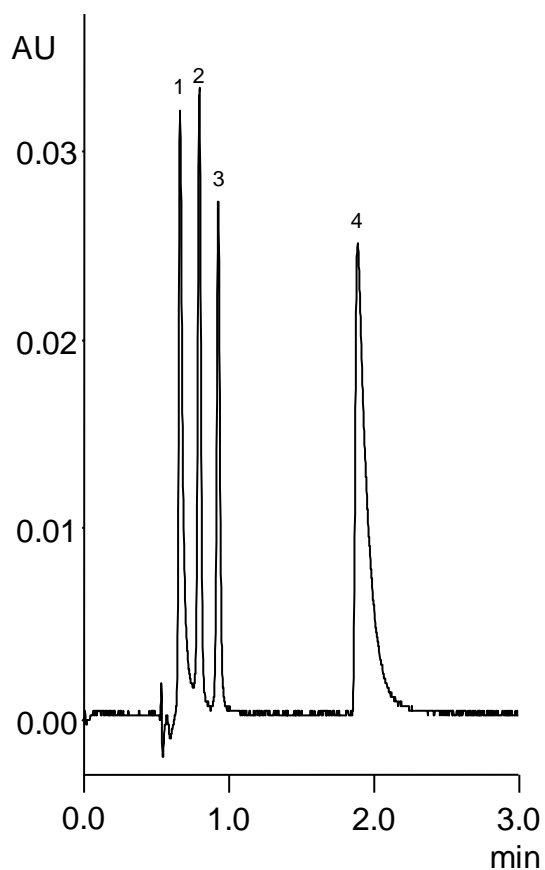


1.

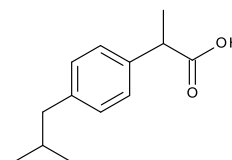


Ibuprofen

2.

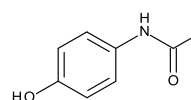


Acetaminophen

3.

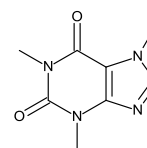

Caffeine

4.

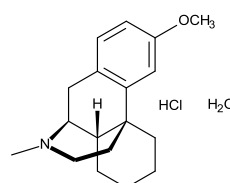

Dextromethorphan HCl H₂O

Column	: Alcyon SFC Triart PFP (5 μm, 12 nm) 150 X 4.6 mm I.D.
Eluent	: CO ₂ /methanol (0.1%DEA) (80/20)
Flow rate	: 3.0 mL/min
Temperature	: 40°C
Detection	: UV at 254 nm
Back Pressure	: 13.8 MPa (2000 psi)
Injection	: 1 μL (0.044mg/mL ~ 5.32 mg/mL)