



レムデシビル
Remdesivir

E210219F

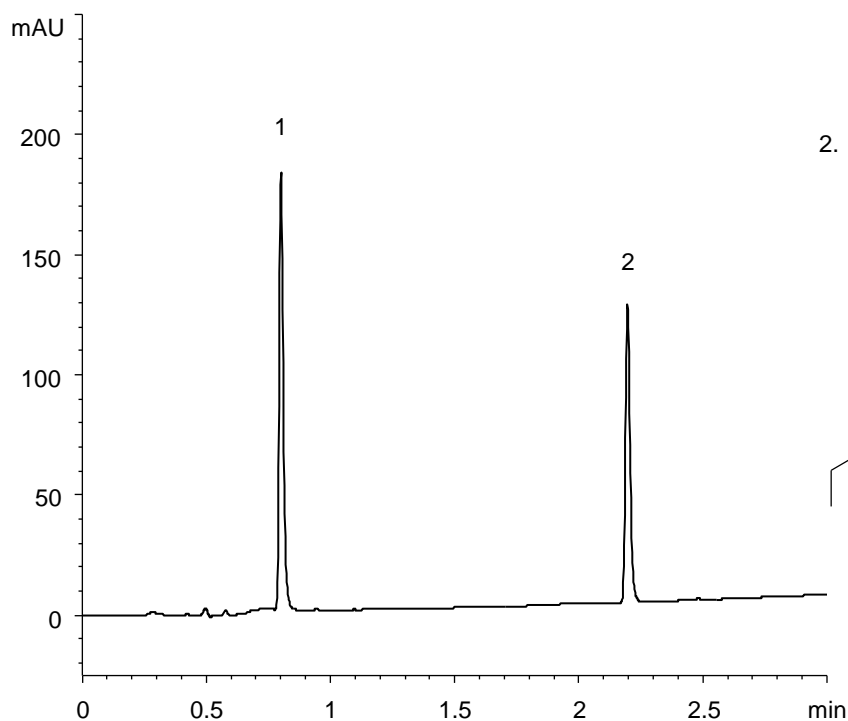


1.

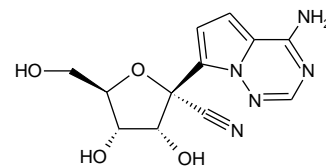

GS-441524

2.

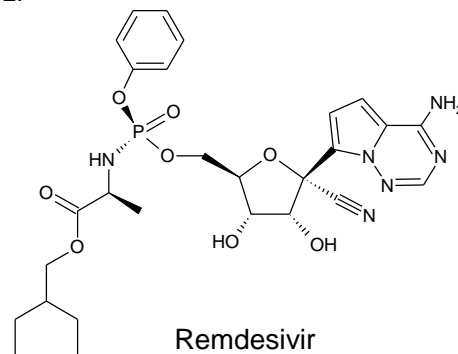

Remdesivir

Column	: YMC-Triart C18 (1.9 μ m, 12 nm) 50 X 2.1 mm I.D.
Eluent	: A) water/formic acid (100/0.1) B) acetonitrile/formic acid (100/0.1) 5-90%B (0-3 min)
Flow rate	: 0.4 mL/min
Temperature	: 40°C
Detection	: UV at 254 nm
Injection	: 2 μ L (0.025 mg/mL, 0.05 mg/mL)