



環状ペプチド
Cyclic peptides


G210901A

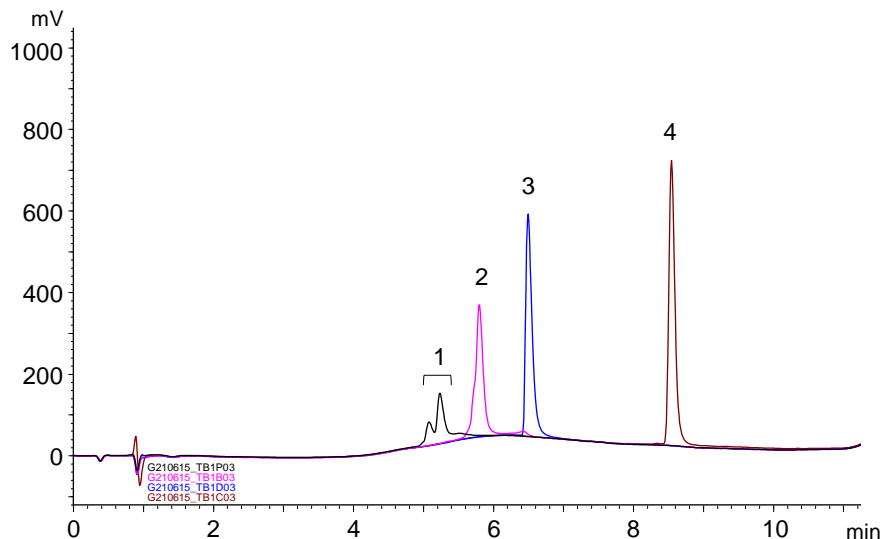


1. Polymyxin B sulfate

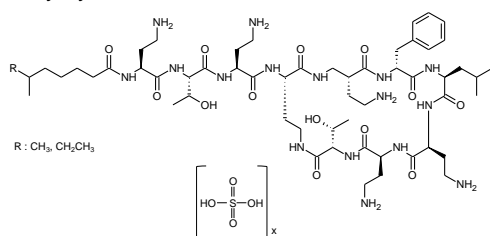


3. Daptomycin

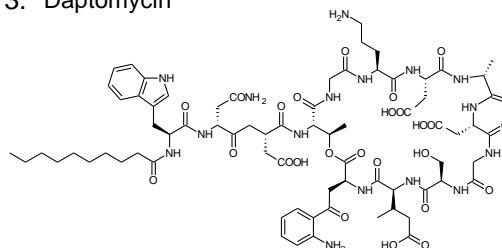

2. Bacitracin

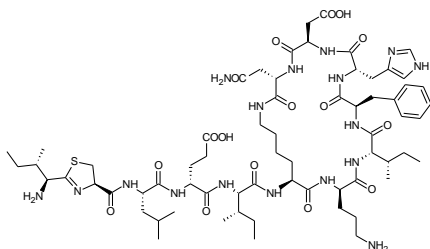

4. Cyclosporin A

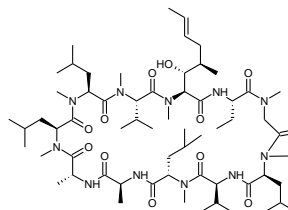

Column	: YMC-Triart Bio C4 (1.9 μm, 30 nm) 50 X 2.1 mm I.D.
Eluent	: A) water/formic acid (100/0.1) B) acetonitrile/formic acid (100/0.1) 5%B (0-0.5 min), 5-90%B (0.5-7.5 min), 90%B (7.5-9 min)
Flow rate	: 0.2 mL/min
Temperature	: 70°C
Detection	: UV at 210 nm
Injection	: 1 μL (0.5 mg/mL)